

UROSKOP D

SP

Log book

Installation of service software

© Siemens AG 1998

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights created by patent grant or registration of a utility model or design, are reserved.

Register 10

Print number: RLL5-310.113.01.05.02

Replaces: RLL5-310.113.01.04.02

English

Doc. gen. date: 11.98

NOTICE

Due to the technologies used in the UROSKOP D, only service engineers who have completed training course LL 04.01.00 may perform service work on the system.

Required documents

- Service PC manufacturer's handbook.
- Operating the service software: RLL5-310.113.02...

Required auxiliary devices

- Service PC
- Service software
- Host software
- Memory software (not for UROSKOP D3)

NOTICE

Service software VA00A and VA00B differ in the file transfer only. Therefore, the instructions for the service software RLL5-310.113.02... apply to both versions.
Service software VA00A will not function on fast (Pentium) PC's.

General

The software is supplied on three diskettes. Two of the diskettes contain the system software for operating the UROSKOP D. The third diskette contains the actual service software:

1. Diskette with service software
INITDISK.EXE
SERVICE.EXE
SERVICE.HLP
README.TXT
System-specific files (*.Cfg, *.Par, etc.).
2. Diskette with Host software
UROS.HEX
3. Diskette with MEMOSKOP 100 software (UROSKOP D1/ D2 only)
BILD.HEX
HARD.HEX
4. Diskette with MEMOSKOP 2K software (UROSKOP D1 only) with Polydoros SU
BILD.HEX

Installing the software under DOS

Close any programs running, until you receive the system message:

```
C:\ >
```

Now, create a directory in which the software can be saved (e. g. UROSKOP D).
Confirm this with the command:

```
C:\ > md UROSKOPD      < ENTER >
```

Select the directory:

```
C:\ > cd UROSKOPD      < ENTER >
```

Subsequently, the message appears:

```
C:\UROSKOPD >         < ENTER >
```

Now insert the diskettes sequentially in the PC internal drive.
This should be named drive A. Select the drive with the following command:

```
C:\UROSKOPD > a:       < ENTER >
```

Now, the message appears:

```
A:\ >
```

The files from the diskette will be copied to the hard disk. In addition, enter:

```
A:\> copy a:\*.* c:\UROSKOPD < ENTER >
```

All files have now been copied from the service software, host software and MEMOSKOP software diskettes to the hard disk.

Now the service software can be started. To do this, use the following command:

```
A:\> c:                < ENTER >
```

Now the system will switch back to drive C. The message appears:

```
C:\ > UROSKOPD
```

After entering:

```
C:\UROSOPD > service < ENTER >
```

the menu for the service software appears.

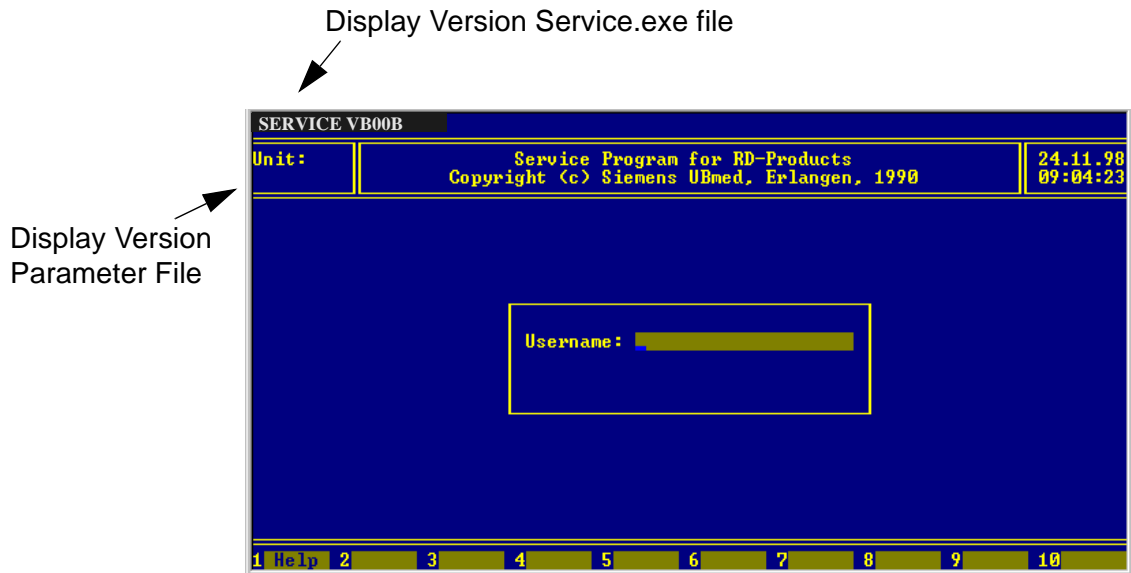


Fig. 1

To exit the program at this point, press the < Ctrl > and < C > keys simultaneously.

Installing the software under Windows 95/ 98

- Open Explorer.
- Select drive "c".
- Click on drive "c" with the right mouse key.
- Click on New -> Folder.
- Enter the folder name (UROSOPD).
- Insert the diskette with the service software in drive "a".
- Select drive "a".
- Highlight all the files.
- Click on Edit -> copy.
- Double-click on drive "c" -> Folder UROSOPD.
- Click in the UROSOP D folder with the right mouse key.
- Select "Paste".
 - The files will be copied from the diskette to the UROSOP D folder.
- Exit Explorer.

Creating an icon on the Desktop

The UROSKOP service software will run under Windows '95 and Windows '98 only in MS-DOS mode.

To create an icon, click with the right mouse key on the desktop. In the menu next to "New", select "Shortcut". Either click on "Browse" and select the service.exe file in the directory (UROSKOPD) or enter the path (c:\uroskopd\service.exe) in the command line. The options /c for color or /s for simulation or, if needed, both, can also be applied. It is advisable to create a separate icon for the simulation. After clicking on "continue", you can enter a name for the shortcut (e.g. UROSKOP service). Select an icon in the next field and click on "Finish".

Click on the icon just created with the right mouse key and select "Properties". Under "Program", select "Advanced". The window "Advanced program settings" opens. In this window, select "MS-DOS mode" (a check mark appears in the box). This ensures that only the service program will run and that all Windows programs will be ended. Exit the "Advanced Program Settings" dialog box with OK and exit the "Properties" window with OK.

Operation under Windows NT and Windows CE

The UROSKOP service software cannot be operated on Windows NT or Windows CE!

TD PS 24 / Erlwein
TD SP / M. Heinlein
SMS Iselin/ O'Donnell

This page intentionally left blank.